BabyOno®

PL

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catalogue no. 291

ENG

The electronic baby scale is intended for monitoring the baby's weight. It is especially recommended in the first months of life and for premature babies whose weight is a particularly important indicator of their growth.

The special design of the scale structure makes it possible to weigh both babies (up to 20 kg) and older children (up to 50 kg), which enables you to use the scale for years. We took care to ensure the highest quality of the product, which guarantees comfort and safety during the weighing.

With your needs in mind, we created a product which is easy to use, accurate and reliable. The scale is equipped with the memory function which makes it possible to compare the current weight of the child with the previous weighing result, as well as the tare weight function available for the whole measuring range. The scale automatically freezes the weighing result, which makes it possible to read it without haste.

The scale is suitable for weighing:

- Babies (maximum load 20 kg)
- Children (maximum load 50 kg)

TECHNICAL SPECIFICATION

Maximum load – 50 kg (the scale for children) and 20 kg (the tray for babies)

Minimum load – 0.05 kg

Gradation – 0.01 kg

Units of measurement – kg/ lb-oz

Automatic turn-off

Automatic freezing of the measurement result

Tare weight function

Memory of the last measurement result

Low battery indicator

Power supply – 2 x AA (1,5 V)


OPERATION MANUAL

Before the first use of the device please make sure to read the operation manual carefully and observe the warnings contained therein. Please pay particular attention to the safety recommendations, so as to avoid any accidents and/or product damage when using the device. Please retain this manual for future reference!

SAFETY MEASURES

1. Use the scale in accordance with its intended use.
2. Do not overload the scale. The overloading is signalled by displaying „O-Ld“ on the LCD screen.
3. Maximum load:
 - 20 kg on the baby tray
 - 50 kg on the scale pan for children
4. Do not expose the scale to shock or impact.
5. Do not immerse the scale in water.
6. Do not place other items on the scale during storage.
7. Do not use or store the scale in a damp environment.
8. The optimal temperature of the environment for the scale operation and storage: from 5°C to 30°C

DESIGN OF THE SCALE

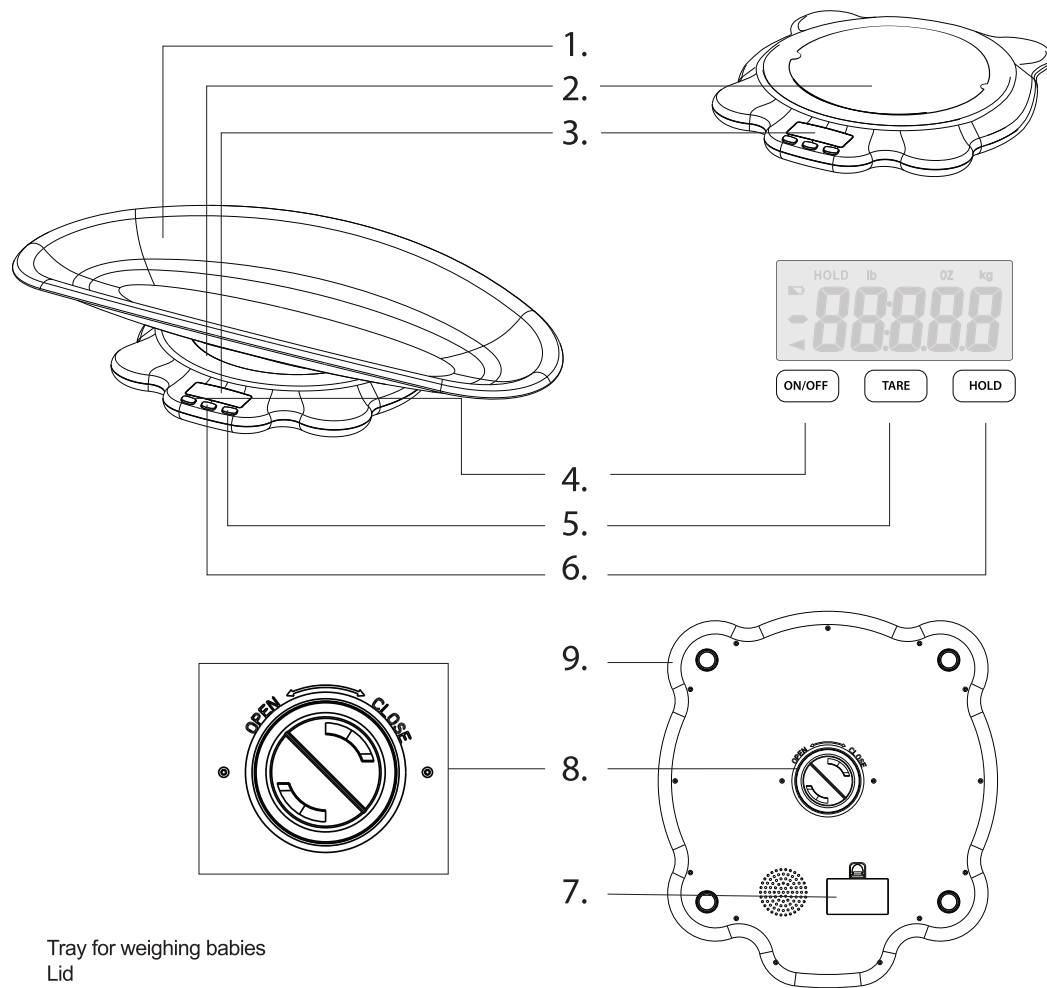Tray for weighing babies

Lid

Display

Power button

Tare button

Button to unfreeze the reading

Battery cover

Safety nut to protect the scale pans: the baby tray (1) and the lid (2)

Rubber feet

SCALE PAN ASSEMBLY INSTRUCTION

1. Unlock the safety nut (8) by turning it to the left to the "Open" position.
2. Place the appropriate scale pan in the recess (depending on the intended use and the maximum load):
 - a) tray (1) for weighing babies and children in lying position – maximum load up to 20 kg
 - b) lid (2) for weighing children in standing position – maximum load up to 50 kg
3. Lock the scale pan with the safety nut (8) by turning it to the right to the "Close" position.
4. Check whether the elements are connected in a correct and stable way.

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OPERATION MANUAL

1. Make sure the scale is placed on a horizontal, stable, dry and hard surface. Do not place the scale on a carpet or flexible surfaces.
2. Press the power button (4) to switch the power supply on.
3. If you press the power button (4) once again before “0.00” appears on the display, you will be able to change the unit of measurement: kg or lb.
4. If you want to set the tare weight, place the object(s) whose weight has to be deducted from the child's weight (e.g. a towel or dressing gown) and press the tare button (5) to zero the scale out.
Note: Since the measuring range starts from 50 g, the scale will not account for the tare weight below this value (e.g. in the case of a thin baby garment).
5. The scale is ready to use when the display shows “0.00”.
6. Place on the appropriate scale pan:
 - a) a baby in lying position, if the baby tray (1) has been mounted – maximum load of 20 kg
 - b) a child in standing position if the scale lid (2) has been mounted – maximum load of 50 kg
7. Wait a few seconds for the scale to establish and display the weighing result. The reading will freeze automatically.
8. When the display shows the symbol “HOLD”, which means the weighing result has been frozen, take the child off the scale.
9. The “HOLD” button unfreezes the reading, which makes it possible to use the scale for another weight measurement. The next weighing may be done when „0.00” appears on the display.
10. Turn off the device by pressing the power button. The device will also turn off automatically after 30 seconds of idleness.
11. The scale has a built-in memory holding the last measurement. When the scale is switched on again, the display will show the last weighing result (Note: the value of the last weighing result should be greater than 0,5 kg).

BATTERY REPLACEMENT

When scale display shows the symbol „” it signals that the battery is low. In order to ensure correct operation of the device, please replace the old batteries with new ones. The battery compartment is located in the bottom of the scale unit.

In order to replace the batteries:

1. Open the battery cover (7)
2. Insert two 1.5V AA batteries into the compartment, following the polarity markings shown inside.
3. Close the battery compartment by pressing the cover.
4. Used batteries should be disposed of in accordance with the currently binding legal regulations.

CLEANING AND MAINTENANCE

Wipe the scale with a damp cloth.
Do not immerse the scale in water!
Do not use any corrosive chemical agents to clean the scale.

The device should be disposed of in accordance with the currently binding legal regulations.

The device is compliant with the requirements of the EU directive: Electromagnetic Compatibility (EMC) 2004/108/CE.

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